

T5 Tricycles



Tracker T5 24"



TMX T5

Features

- Carefully engineered frames for strength and style
- Low step-through frame for easy access
- A choice of 3 frame sizes
- Front suspension fork for improved performance and handling
- Front Disc brake with additional V-style brake with parking option for ultimate stopping power
- 3 gearing options – Single speed fixed-wheel, single speed free-wheel & 5 Speed Shimano indexed gear
- Adjustable steering stop to assist stability
- Many optional accessories which allow the tailoring of the tricycle to suit every rider's needs

A whole new world of possibilities



Derived from Theraplay's popular range of foot-propelled tricycles, the T5 Tricycle is designed for the more adventurous cyclist. Created with a robust frame that can absorb the bumps and hollows of our streets and pathways, the T5 gives the rider the ability to explore wider horizons with greater confidence.


Theraplay's design team have given special attention to the styling of the T5 Tricycle. Care has been taken to select components that will give maximum performance and enhance the appearance of the tricycle. Designed to suit the needs of teenage to adult riders, the T5 model incorporates a comfortable seating position that adjusts for growth, variable height handlebars and two forms of positive easy to use brakes.

T5 Tricycles are available with a choice of three types of gearing. A 5 Speed Shimano Indexed gearing system gives the potential for longer distance rides where the rider is able to select the gear ratio depending upon the desired speed or gradient of hill. A single speed free-wheel gear is for riders who wish to cycle without the pedals constantly turning. This option will allow riders to navigate cycle paths and parks with ease. A fixed-gear provides the maximum therapeutic benefit for the rider and gives increased manoeuvrability. The foot pedals turn as the wheels turn, meaning both legs are given a thorough workout. This type of gearing enables the rider to control the speed of the tricycle with the foot pedals and to reverse if caught in a tight corner.

Optional Accessories

- Frame Folding Mechanism
- Adjustable Back/Trunk Support
- Foot Support Systems
- Saddle Options
- Loop Style or Cumfigrip Handlebars
- Pneumatic or Puncture Proof Tyres
- Rear Carrier Basket
- Carer Control Pole

Colours

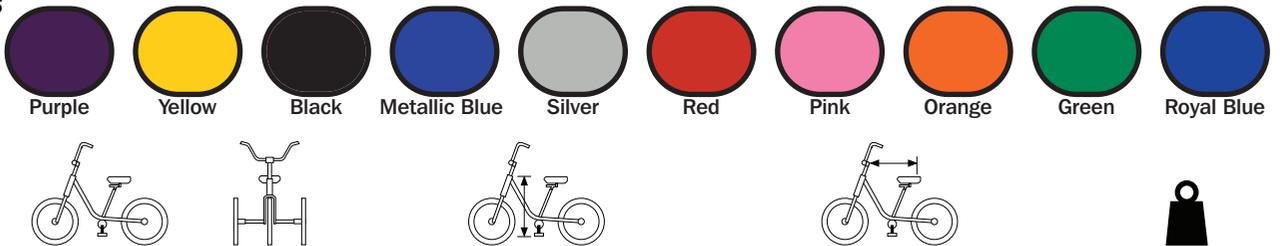

	Overall Length	Overall Width	To suit inside leg		Seat to handlebars		Max user weight
			(MIN)	(MAX)	(MIN)	(MAX)	
TMX	59" (1500mm)	32" (810mm)	20" (510mm)	25" (640mm)	19" (480mm)	22" (560mm)	130lbs (64kg)
Tracker 20"	65" (1650mm)	32" (810mm)	24" (610mm)	29" (740mm)	21" (535mm)	24" (610mm)	170lbs (77kg)
Tracker 24"	69" (1753mm)	32" (810mm)	25" (640mm)	30" (762mm)	22" (560mm)	25" (640mm)	170lbs (77kg)



32 Welbeck Road, Darnley Industrial Estate, Glasgow, G53 7SD, Scotland
 Tel: 0141 876 9177 Fax: 0141 876 9039
www.theraplay.co.uk



PO Box 1364, Cumberland, MD 21501-1364, USA
 Tel: 301 759 3525 or 800 306 6777 Fax: 301 759 3525
www.triaid.com

Local Dealer

